Levenhuk Kelvin TMC13 LED Microscope for Smartphone

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
levenhuk
Zoom&Joy

OPTICS
FOR
SMARTPHONES **Kelvin**
MASTER THE KEY LIGHT

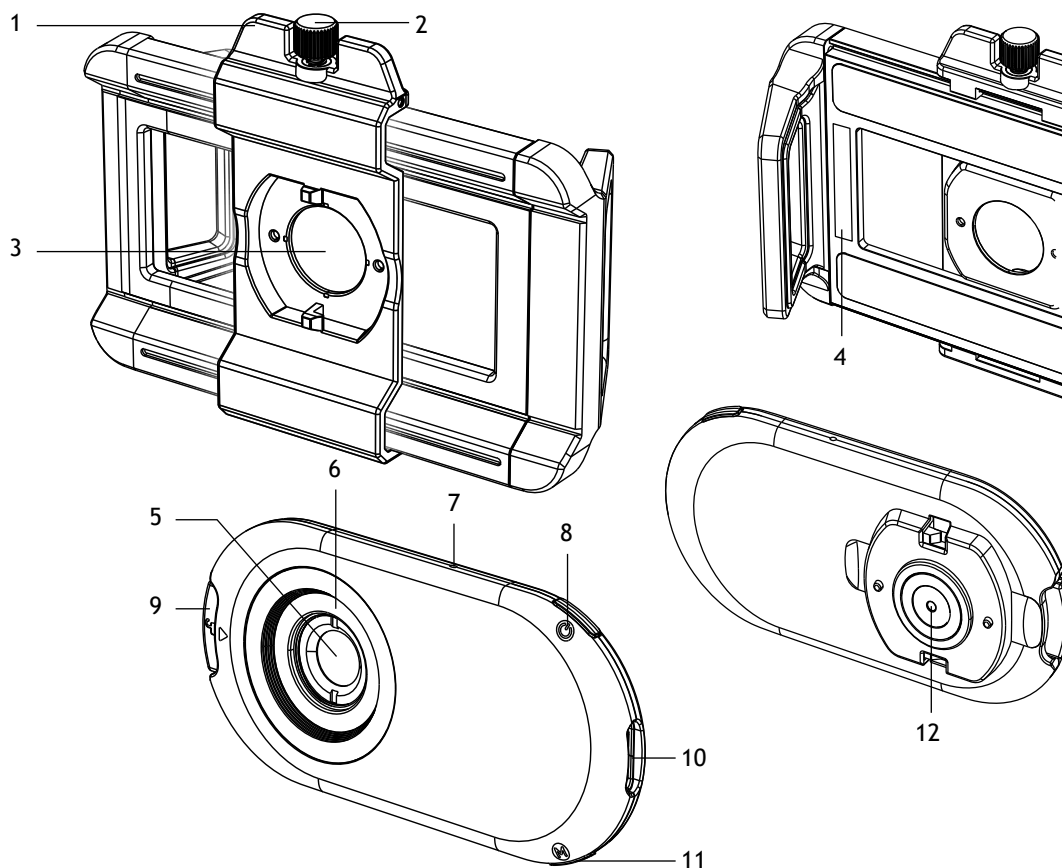EN	BG	CZ	DE	ES
1 Movable lens mount	Подвижна монтировка за обектив	Posuvné uchycení objektivu	Verstellbare Objektivhalterung	Montura de lente móvil
2 Fixing screw	Застопоряващ винт	Upevňovací šroub	Befestigungsschraube	Tornillo de fijación
3 Lens hole	Отвор за обектива	Otvor pro objektiv	Objektivöffnung	Orificio para la lente
4 Placement area for EVA pads	Място за поставяне на подложки от EVA	Místo pro umístění podložek EVA	Auflagefläche für EVA-Polster	Área de colocación de los espaciadores de EVA
5 Objective lens	Леща на обектива	Čočka objektivu	Objektivlinse	Lente objetivo
6 LED illumination	Светодиодно осветление	LED osvětlení	LED-Beleuchtung	Iluminación LED
7 Power indicator	Индикация за захранването	Indikace napájení	Leistungsanzeige	Indicación de encendido
8 Power button	Захранващ бутон	Tlačítko napájení	Ein-/Aus-Taste	Botón de encendido
9 Unlock button	Бутон за отключване	Tlačítko odemknutí	Entriegelungstaste	Botón de desbloqueo
10 USB Type-C slot	Слот за USB Тип-C	Slot pro USB typu C	USB-Typ-C-Buchse	Puerto USB de tipo C
11 LED illumination mode button	Бутон за режим на светодиодно осветление	Tlačítko režimu LED osvětlení	LED-Beleuchtungsmodus-Taste	Botón del modo de iluminación LED
12 Eyepiece	Окуляр	Okulár	Okular	Ocular

FR	HU	IT	PL	PT	RU
1 Support d'objectif amovible	Mozgó lencsefoglalat	Montatura mobile per lente	Ruchome mocowanie do obiektywu	Suporte móvel da lente	Подвижное крепление для объектива
2 Vis de fixation	Rögzítőcsavar	Vite di fissaggio	Śruba mocująca	Parafuso de fixação	Крепежный винт
3 Ouverture de l'objectif	Lencse nyílása	Foro lente	Otwór na obiektyw	Abertura da lente	Отверстие для объектива
4 Zone de placement des patins EVA	EVA párnák elhelyezési területe	Area di posizionamento per cuscinetti in EVA	Miejsce na podkładki EVA	Área de colocação de almofadas de proteção EVA	Место установки прокладок из ЭВА
5 Objectif	Objektívlencse	Lente obiettivo	Obiektyw	Lente da objetiva	Объектив
6 Éclairage à DEL	LED-es megvilágítás	Illuminazione LED	Oświetlenie LED	Iluminação LED	Светодиодная подсветка
7 Témoin d'alimentation	Bekapcsolásjelző	Indicatore di accensione	Wskaźnik zasilania	Indicador de energia	Индикатор питания
8 Bouton d'alimentation	Bekapcsológomb	Pulsante di alimentazione	Przycisk zasilania	Botão de ligar/desligar	Кнопка питания
9 Bouton de déverrouillage	Feloldógomb	Pulsante di sblocco	Przycisk odblokowania	Botão de desbloqueio	Кнопка разблокировки
10 Prise USB type C	C-típusú USB-foglalat	Slot per USB Tipo-C	Gniazdo USB typu C	Porta USB tipo C	Разъем USB Type-C
11 Bouton du mode d'éclairage à DEL	LED-es megvilágítási mód gomb	Pulsante modalità illuminazione LED	Przycisk trybu oświetlenia LED	Botão do modo de iluminação LED	Кнопка режима светодиодной подсветки
12 Oculaire	Szemlencse	Oculare	Okular	Ocular	Окуляр

EN Levenhuk Kelvin TMC13 LED Microscope for Smartphone


Caution! Never look directly at the Sun through this device, as this may cause permanent eye damage and even blindness.

The kit includes: LED microscope, phone-to-lens mounting bracket, USB Type-C cable, EVA pads (9 pcs), pouch, cleaning cloth, user manual, and warranty.

Charging the device

Connect the USB Type-C cable to the USB Type-C slot (10) and power bank, PC, or power adapter (not included). The power indicator (7) will light red. Charging takes approximately 90 minutes. When charging is complete, the power indicator will light green.

Getting started

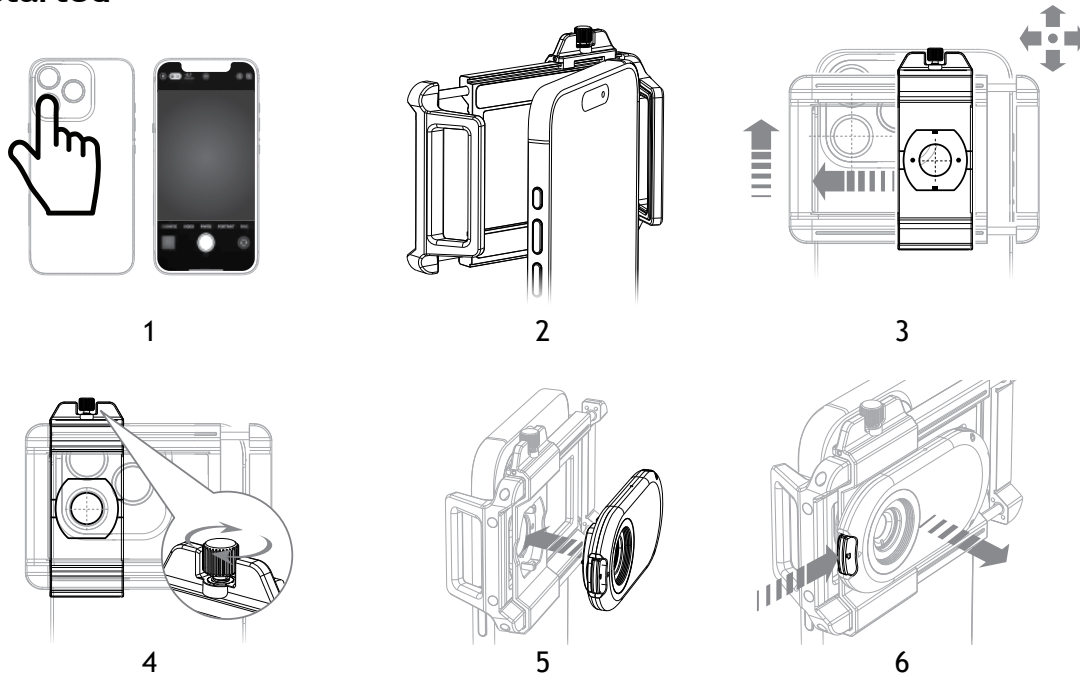- Find the main camera on your smartphone. To do so, open the camera app on your smartphone and cover each lens on the back one by one with your finger. The module that makes the image on the screen go black when covered is the main camera (Fig. 1).
- Gently push your smartphone into the phone-to-lens mounting bracket (Fig. 2). To precisely adjust the mounting bracket to the thickness and width of your phone, use the soft EVA pads (included).
- Move the lens mount (1) horizontally and the mounting bracket vertically to align the lens hole (3) with your smartphone's main camera (Fig. 3).
- Tighten the fixing screw (2) (Fig. 4).
- Install the microscope directly into the lens mount (1) until you hear a click (Fig. 5).
- Press and hold (9) to remove the microscope from the mounting bracket (Fig. 6).

Usage

- Press and hold (8) to turn the device on. The LED illumination (6) will light up.
- Press (8) to change the LED brightness level (low/middle/high).
- Press (11) to change the LED illumination mode. You may choose from 2 color temperatures (cool or warm light) in 3 illumination ranges (full circle / left half circle / right half circle) in each of 3 brightness modes (low/middle/high), or you may select the UV full circle mode.
- Place the microscope directly above the object you want to capture.
- Press and hold (8) to turn the device off.

! If there are uniform dark corners around the image, adjust the digital zoom on your smartphone. If the dark corners look uneven, check that the lens hole (3) is aligned with your smartphone's main camera before trying again.

Specifications

Magnification	300x
Focusing distance	10mm
Objective lens diameter	22mm
Optics material	optical glass
Optics coating	fully multi-coated
CPL filter	+
Brightness level	3 modes
LED illumination	19 modes
Compatibility with smartphone models	iPhone, Android; smartphone width: 71–105mm (2.8–4.1"), smartphone thickness: ≤ 20mm (0.8")
EVA pads	21x4.5x2mm (3 pcs) 21x4.5x3mm (3 pcs) 21x4.5x4mm (3 pcs)
Operating temperature range	0... +45° C (+32... 113° F)
Power supply	3.7V, 100mA·h rechargeable Li-ion battery; 5V, 1A DC adapter (not included); USB Type-C cable (included)
Battery duration	2 h (1 h with UV lights on)

The manufacturer reserves the right to make changes to the product range and specifications without prior notice.

Care and maintenance

Take the necessary precautions when using the device with children or others who have not read or who do not fully understand these instructions. Do not try to disassemble the device on your own for any reason. Protect the device from water, dust, sudden impact and excessive mechanical force. Do not touch the optical surfaces with your fingers. Clean the lens surface with compressed air or a soft lens cleaning cloth. Never attempt to operate a damaged device or a device with damaged electrical parts! Store the device in a dry cool place. **If a part of the device or battery is swallowed, seek medical attention immediately.**

Battery safety instructions

The device is equipped with a rechargeable lithium-ion battery. This avoids frequent battery replacement. Always switch the device off when not in use. If the battery charge is low, please recharge the device in time. Do not overheat the battery. Do not discharge the battery completely. Do not disassemble, puncture, or deform the battery. Do not expose the battery to fire, high temperatures, or immerse it in water or other liquids. Improper handling of the battery may result in fire, explosion, or burns. If the battery becomes damaged, deformed, or swollen, discontinue use immediately. Keep batteries out of the reach of children, to avoid risk of ingestion, suffocation, or poisoning. Utilize used batteries as prescribed by your country's laws.

Levenhuk International Warranty

All Levenhuk telescopes, microscopes, binoculars and other optical products, except for accessories, carry a **50-year warranty** against defects in materials and workmanship. This warranty is valid for **50 years** from the date of purchase. All Levenhuk accessories are warranted to be free of defects in materials and workmanship for six months from date of retail purchase. The warranty entitles you to free repair or replacement of the Levenhuk product in any country where a Levenhuk office is located if all warranty conditions are met.

For further details, please visit: levenhuk.com/warranty

If warranty problems arise, or if you need assistance in using your product, contact the local Levenhuk branch.

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Levenhuk® is registered trademark of Levenhuk, Inc.
© 2006–2026 Levenhuk, Inc. All rights reserved.
Levenhuk Inc. (USA): 6021 Catlin Dr., Tampa, FL 33612, USA,