
Levenhuk Kelvin HPM10
Kids Photomicroscope


EN User Manual
BG Ръководство за потребителя
CZ Návod k použití
DE Bedienungsanleitung
ES Guía del usuario
FR Mode d’emploi
HU Használati útmutató
IT Guida all’utilizzo
PL Instrukcja obsługi
PT Manual do usuário
RU Инструкция по эксплуатации

levenhuk
Zoom&Joy

OPTICS FOR SMARTPHONES Kelvin
MASTER THE KEY LIGHT

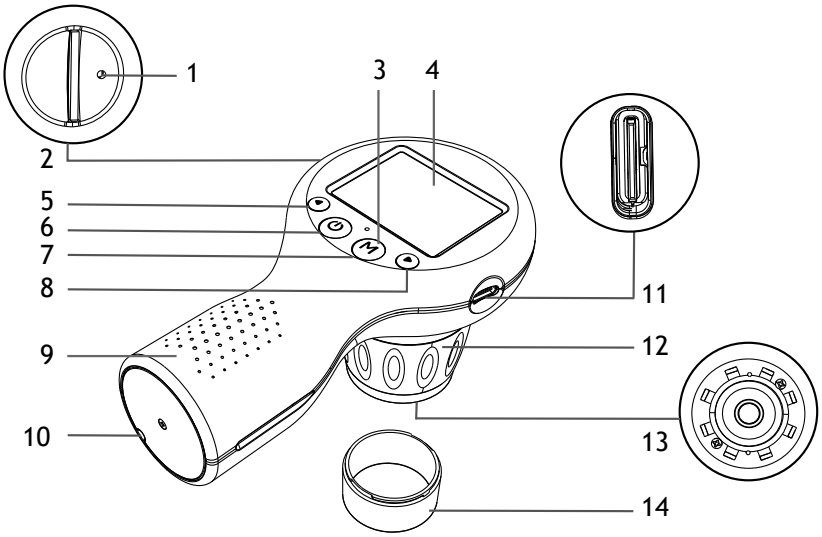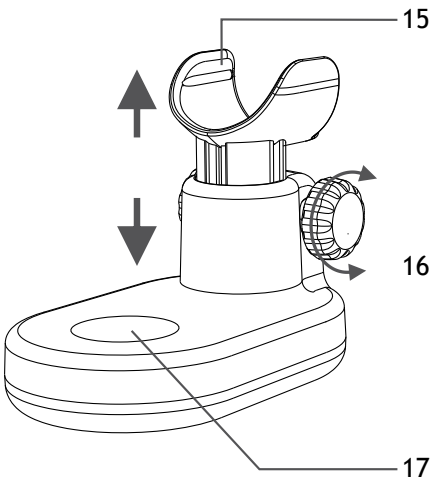	EN	BG	CZ	DE	ES
1	Reset button	Бутон Възстановяване	Tlačítko Resetovat	Wiederherstellen-Taste	Botón Restablecer
2	MicroSD card slot	Слот за microSD карта	Slot pro kartu microSD	MicroSD-Kartensteckplatz	Ranura para tarjeta microSD
3	Power indicator	Индикация за захранването	Indikace napájení	Leistungsanzeige	Indicación de encendido
4	Screen	Екран	Obrazovka	Display	Pantalla
5	◀ (Brightness down/ Zoom out/Left) button	Бутон ◀ (Яркост надолу/ Отдалечаване/Наляво)	Tlačítko ◀ (Snížení jasu/ Oddálení/Vlevo)	◀ -Taste (Helligkeit verringern/ Vergrößerung aus-/Links)	Botón ◀ (Disminuir brillo/ Alejar/Izquierda)
6	⏻ (Power/Confirm/Photo) button	Бутон ⏻ (Захранване/ Потвърждаване/Снимка)	Tlačítko ⏻ (Napájení/ Potvrzení/Foto)	⏻ -Taste (Ein-/Aus-/Bestätigen/ Foto)	Botón ⏻ (Encendido/ Confirmar/Foto)
7	Ⓜ (Mode/Menu) button	Бутон Ⓜ (Режим/Меню)	Tlačítko Ⓜ (Režim/Nabídka)	Ⓜ -Taste (Modus/Menü)	Botón Ⓜ (Modo/Menú)
8	▶ (Brightness up/ Zoom in/Right) button	Бутон ▶ (Яркост нагоре/ Приближаване/Дясно)	Tlačítko ▶ (Zvýšení jasu/ Přiblížení/Doprava)	▶ -Taste (Helligkeit steigern/ Vergrößerung ein-/Rechts)	Botón ▶ (Aumentar brillo/ Acercar/Derecha)
9	Handle	Ръкохватка	Rukojeť	Griff	Mango
10	Lanyard hole	Отвор за ремък	Otvor pro šňůrku	Lanyard-Loch	Orificio para cordón
11	USB Type-C slot	USB слот Тип-C	Slot USB typu C	USB-Typ-C-Anschluss	Puerto USB de tipo C
12	Focusing ring	Пръстен за фокусиране	Zaostřovací kroužek	Fokussierring	Anillo de enfoque
13	Objective lens	Обектив	Objektiv	Objektiv	Objetivo
14	Transparent plastic tube	Прозрачна пластмасова тръба	Průhledná plastová trubice	Transparenter Kunststoffzylinder	Tubo de plástico transparente
15	Microscope holder	Държач за микроскоп	Držák mikroskopu	Mikroskophalter	Soporte del microscopio
16	Height adjustment knob	Копче за регулиране на височината	Knoflík pro nastavení výšky	Höhenverstellknopf	Perilla de ajuste de altura
17	Specimen stage	Предметна маса	Stolek pro vzorky	Objekttisch	Platina para muestras

	FR	HU	IT	PL	PT	RU
1	Bouton de réinitialisation	Alaphelyzetbe állítás gomb	Pulsante Reimposta	Przycisk resetowania	Botão Repor	Кнопка сброса
2	Logement de la carte microSD	Micro-SD-kártyafoglalat	Slot per microSD	Gniazdo karty pamięci microSD	Ranhura de cartão microSD	Разъем для карты microSD
3	Témoin d'alimentation	Bekapcsolásjelző	Indicatore di accensione	Wskaźnik zasilania	Indicador de energia	Индикатор питания
4	Écran	Kijelző	Schermo	Wyświetlacz	Ecrã	Экран
5	Bouton ◀ (Réduction de la luminosité / Zoom arrière / Gauche)	◀ (Fényerő csökkentése/ Kicsinyítés/Bal) gomb	Pulsante ◀ (Diminuzione luminosità/Zoom indietro/ Sinistra)	Przycisk ◀ (Zmniejszenie jasności/Oddalenie/ W lewo)	Botão ◀ (Reduzir luminosidade/Diminuir zoom/Esquerda)	Кнопка ◀ (Уменьшение яркости/Отдаление/ Влево)
6	Bouton ⏻ (Alimentation / Confirmation / Photo)	⏻ (Főkapcsoló/ Megerősítés/ Fényképezés) gomb	Pulsante ⏻ (Alimentazione/ Conferma/Foto)	Przycisk ⏻ (Zasilanie/ Potwierdzenie/Zdjęcie)	Botão ⏻ (Ligar/ Desligar/Confirmar/ Fotografia)	Кнопка ⏻ (Питание/ Подтверждение/Фото)
7	Bouton Ⓜ (Mode / Menu)	Ⓜ (Mód/Menü) gomb	Pulsante Ⓜ (Modalità/ Menu)	Przycisk Ⓜ (Tryb/Menu)	Botão Ⓜ (Modo/Menu)	Кнопка Ⓜ (Выбор режима / Меню)
8	Bouton ▶ (Augmentation de la luminosité / Zoom avant / Droite)	▶ (Fényerő növelése/Nagyítás/ Jobb) gomb	Pulsante ▶ (Aumento luminosità/Zoom avanti/ Destra)	Przycisk ▶ (Zwiększenie jasności/Przybliżenie/ W prawo)	Botão ▶ (Aumentar luminosidade/Aumentar zoom/Direita)	Кнопка ▶ (Увеличение яркости/Приближение/ Вправо)
9	Poignée	Fogantyú	Impugnatura	Uchwyt	Pega	Ручка
10	Trou de sangle	Zsinór lyuk	Foro per cordino da collo	Otwór na smycz	Orifício para cordão de segurança	Отверстие для ремешка
11	Prise USB type C	C-típusú USB foglalat	Slot per USB Tipo-C	Gniazdo USB typu C	Porta USB tipo C	Разъем USB Type-C
12	Bague de mise au point	Fókuszáló gyűrű	Ghiera della messa a fuoco	Pierścień regulacji ostrości	Anel de focagem	Колесо фокусировки
13	Objectif	Objektív	Obiettivo	Obiektyw	Objetiva	Объектив
14	Tuyau plastique transparent	Átlátszó műanyag cső	Tubo di plastica trasparente	Przerzoczysta rurka plastikowa	Tubo de plástico transparente	Прозрачный пластиковый тубус
15	Support du microscope	Mikroszkóptartó	Portamicroscopio	Uchwyt mikroskopu	Suporte do microscópio	Держатель микроскопа
16	Bouton de réglage de hauteur	Magasságállító gomb	Manopola di regolazione altezza	Pokrętło regulacji wysokości	Botão de regulação de altura	Ручка регулировки высоты
17	Platine d'observation	Tárgyasztal	Tavolino dei campioni	Stolik na preparaty	Platina de espécime	Предметный столик

EN Levenhuk Kelvin HPM10 Kids Photomicroscope

Caution! Never look directly at the Sun through this device, as this may cause permanent eye damage and even blindness.

The kit includes: digital photomicroscope, base, USB Type-C cable, lanyard, specimen container, cleaning cloth, user manual, and warranty.

Charging the device

This device uses a rechargeable Li-ion battery. Connect the USB Type-C cable to the USB Type-C slot (9) and power bank, PC, or power adapter (not included). The power indicator (3) will light red. Charging takes approximately 1.5 hours. When charging is complete, the power indicator will light green.

! Charging must be performed under adult supervision only. Use only a 5V/1A power source.

Getting started

- Press and hold the ⏻ button (6) to turn the device on. The screen (4) will light up.
- Press and hold the ⏻ button (6) to turn the device off.

Installing/removing the memory card

The device supports microSD memory cards up to 128GB (not included).

- Insert the memory card (not included) into the memory card slot (2). Make sure the direction is right. Do not force the card in.
- Format the memory card before starting, as it will improve its compatibility with the device.
- To remove the card from the slot, press the edge of the memory card in gently, and the card will pop out. Then, remove the card.

Mode selection

Press the Ⓜ button (7) to switch between modes (photo, video, and playback). Photo mode is the default setting.

Settings

Press and hold the Ⓜ button (7) to enter the settings, and then press the ▶ button (8) to choose the settings menu (photo, video, general). Press the ⏻ button (6) to view the settings. Use the ◀ button (5) to scroll up and the ▶ button (8) to scroll down, then press the ⏻ button (6) to confirm.

You may now use the device to observe objects, take photos, and record videos.

Usage

- Aim the device at the object to be observed or ensure that the object is correctly positioned inside the transparent plastic tube (14) below the objective lens (13) (Fig. 1).
- Rotate the focusing ring (12) until the view is in focus (Fig. 2).

Base

- Place the microscope handle (9) into the holder (15) (Fig. 3).
- Place the object to be observed on a specimen stage (17). Rotate the height adjustment knob (16) until the view is in focus (Fig. 4).

Photo and video

Press the ⏻ button (6) to take a photo or to start recording a video. Press the ▶ button (8) to zoom in and the ◀ button (5) to zoom out. Press and hold the ▶ button (8) to turn the LED illumination on and increase brightness, and the ◀ button (5) to decrease brightness and turn the LED illumination off.

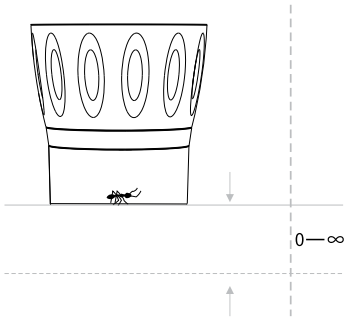
Fig. 1

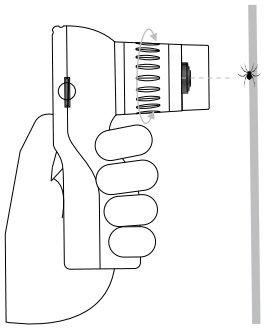
Fig. 2

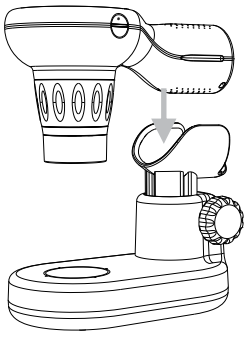
Fig. 3

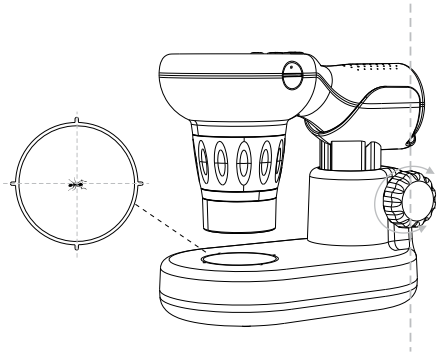
Fig. 4

Playback
In playback mode, press the ◀ button (5) or the ▶ button (8) for scrolling through your pictures and video. Press the ⏻ button (6) to playback a video file. Press and hold the ◀ button (5) to delete the current file.

Copying files
You can transfer files from the device to your PC directly.

- Insert the memory card into a card reader; or
- Connect the device to your computer with a USB Type-C cable.

Software
To work with the device on a PC, you need to install the VideoCap software, available for download on the official Levenhuk website. This software enables you to view your photos and videos in real time on your PC screen.

Specifications

Research method	bright field
Magnification	microscope: 3.3–20x objective lens: 0–20x
Min. size of the objects to be seen	3μm
Field of view	2.3–7.8mm
Optics material	optical glass
Optics coating	multi-coated
Objective lens type	achromatic
Microscope head	digital screen / PC screen
Display	2" color IPS screen, 320x240px
Sensor	CMOS, 720p
Photo format and resolution	*.jpg: 1MP, 2MP
Video format and resolution	*.avi: 720p
Brightness level	3 modes (high/low/off)
LED illumination	8 pcs white LEDs
Display auto-off	30s, 1 min, 2 min
Device auto-off	3 min, 5 min, 10 min
Storage	microSD card up to 128GB (not included)
Software requirements	Windows 7 and above
Operating temperature range	0… +45°C (+32… 113°F)
Power supply	3.7V, 800mA·h rechargeable Li-ion battery; 5V, 1A DC adapter (not included); USB Type-C cable (included)
Battery duration	2 h

The manufacturer reserves the right to make changes to the product range and specifications without prior notice.

Care and maintenance
Take the necessary precautions when using the device with children or others who have not read or who do not fully understand these instructions. Do not try to disassemble the device on your own for any reason. Protect the device from water, dust, sudden impact and excessive mechanical force. Do not touch the optical surfaces with your fingers. Clean the lens surface with compressed air or a soft lens cleaning cloth. Never attempt to operate a damaged device or a device with damaged electrical parts! Store the device in a dry cool place. **If a part of the device or battery is swallowed, seek medical attention immediately.**

Battery safety instructions
The device is equipped with a rechargeable lithium-ion battery. This avoids frequent battery replacement. Always switch the device off when not in use. If the battery charge is low, please recharge the device in time. Do not overheat the battery. Do not discharge the battery completely. Do not disassemble, puncture, or deform the battery. Do not expose the battery to fire, high temperatures, or immerse it in water or other liquids. Improper handling of the battery may result in fire, explosion, or burns. If the battery becomes damaged, deformed, or swollen, discontinue use immediately. Keep batteries out of the reach of children, to avoid risk of ingestion, suffocation, or poisoning. Utilize used batteries as prescribed by your country’s laws.

Levenhuk International Warranty

All Levenhuk telescopes, microscopes, binoculars and other optical products, except for accessories, carry a **50-year warranty** against defects in materials and workmanship. This warranty is valid for 50 years from the date of purchase. All Levenhuk accessories are warranted to be free of defects in materials and workmanship for **six months** from date of retail purchase. The warranty entitles you to free repair or replacement of the Levenhuk product in any country where a Levenhuk office is located if all warranty conditions are met.

For further details please visit our web site: levenhuk.com/warranty

If warranty problems arise, or if you need assistance in using your product, contact the local Levenhuk branch.

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Levenhuk® is registered trademark of Levenhuk, Inc.
© 2006–2026 Levenhuk, Inc. All rights reserved.
Levenhuk Inc. (USA): 6021 Catlin Dr., Tampa, FL 33612, USA,